lavendu

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SONIC KID TOOTHBRUSH
Model: LST8443

TABLE OF CONTENTS

1. INTENDED PURPOSE.....	17
2. TECHNICAL SPECIFICATIONS	17
3. SAFETY INSTRUCTIONS.....	17
4. EXPLANATION OF SYMBOLS.....	19
5. DEVICE COMPONENTS	20
6. SET CONTENTS.....	21
7. OPERATION.....	21
8. CLEANING AND CARE	24
9. REPAIR	25
10. STORAGE AND TRANSPORT	25
11. DISPOSAL.....	25
12. CE DECLARATION.....	26
13. WARRANTY.....	26

1. INTENDED PURPOSE

The sonic toothbrush is designed for home dental care of children 3 years old and above.

It is not intended for industrial or commercial use (e.g. by medical or veterinary establishments).

Please follow these instructions to ensure safe use of the device.

2. TECHNICAL SPECIFICATIONS

Model	LST8443
Power supply	1 x Battery AAA 1,5 V  (included)
IP class	IPX7

3. SAFETY INSTRUCTIONS

1. Read the entire instruction manual before you start using the device.
2. The device is intended for household use only.
3. It is recommended that children not be left unattended when using the toothbrush.
4. Use the device only for its intended purpose.
5. Do not use the device if it is damaged, and do not attempt to repair it.
6. The device is not a toy.
7. Do not store the device in high humidity.

8. Do not store the device in high temperatures or direct sunlight.
9. Protect the device against shock and impact.
10. Never use strong cleaning agents, solvents or similar substances for cleaning the device.
11. Do not use the product on animals.
12. Wash the toothbrush head after each use.
13. Do not exert excessive pressure on teeth or gums while brushing.
14. Do not replace the head with wet hands.
15. Remove the battery before prolonged storage of the device.
16. Keep this manual and, if possible, the package.
17. If the device is passed on to another user, it must be accompanied by this manual.
18. You may request the user's manual in the electronic version at: <http://instrukcje.vershold.com>

4. EXPLANATION OF SYMBOLS

Read this instruction manual.

The product is compliant with the requirements of the applicable EU directives.

Disposal of waste electrical and electronic equipment – see the DISPOSAL section of this manual.

IPX7

The device can be immersed in water for up to 30 minutes at a depth up to 1 m. "Underwater" waterproof quality – first class. Items of this class can be used when swimming in water and during snorkelling, where no great depths occur. The housing may leak when submerged deeper than 1 m, submerged for longer than 30 minutes or where water pressure is very high.


Symbol indicating that the manufacturer has contributed financially to the construction and operation of a packaging material recovery and recycling system.

Designation of the packaging material – corrugated cardboard.

5. DEVICE COMPONENTS

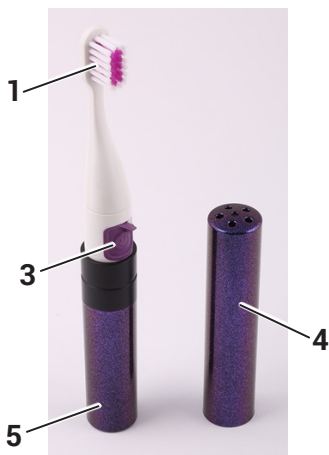1. Head
2. Fastening element
3. On/off switch
4. Top cover
5. battery compartment cover
6. Battery compartment

6. SET CONTENTS

- 1 x Sonic toothbrush
- 1 x Battery AAA 1,5 V --- (included)
- 2 x Head (1 installed, 1 spare)
- 1 x User's manual

Open the packaging and carefully take out the device. Inspect the set for completeness and any damage. Make sure that plastic parts are not cracked. If any part is damaged or missing, contact the dealer and do not use the device. Keep the packaging or dispose of it in accordance with local regulations.


Caution! For the safety of children, do not leave any packaging parts accessible and unattended (plastic bags, cardboard boxes, Styrofoam etc.). Danger of suffocation!

7. OPERATION

7.1 Battery installation/replacement


1. Make sure that the device is turned off and is not working.
2. Turn the battery compartment cover as shown in the figure below.
3. Remove the battery compartment cover.
4. Remove the used battery (if present).
5. Insert a new 1.5 V --- AAA battery into the battery compartment as indicated.
6. Replace the cover of the battery compartment and turn it according to the information in the figure below. sure that the device is off.

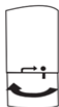Open

Close

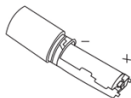Inserting the battery correctly
Overview drawing

Caution!


- The device is designed to be used with AAA 1.5 V --- batteries.
- The battery are included.
- When installing new batteries, observe the markings in the battery compartment.
- Use the type of batteries recommended for this device: AAA 1.5 V --- .
- Do not put used batteries in a rubbish bin, but in special containers for used batteries.
- Never re-charge single-use (non-rechargeable) batteries.
- Do not short-circuit the battery terminals.
- Never combine old batteries with new ones.
- Do not expose batteries to heat (direct sunshine, heaters, fire etc.).
- Risk of explosion! Never disassemble or short-circuit batteries or throw them into fire.
- Always remove the batteries before prolonged storage of the device to prevent damage from possible electrolyte leak.

- Remove dead batteries from the device immediately.
- Dead batteries may leak, which may damage the device.
- If you touch acid from a battery, rinse your hands under running water. If any acid gets into your eyes, contact a doctor. The acid contained in the battery may cause irritation or burns.
- Swallowing a battery may have fatal consequences! Keep batteries away from children and pets. Seek medical assistance immediately if a battery is swallowed.

7.2 Installing / replacing the head

1. Ensure that the device is off.
2. Remove the brush head by turning it clockwise.
3. Replace the new head and turn it over, see picture below.

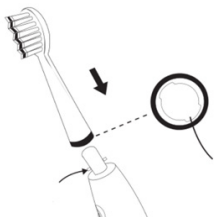The correct application

The correct direction of turning the head during assembly.

Caution!

Warning! The pictorial drawing may differ from the actual product.

Warning!

The head must be clearly latched during installation.

To exclude the risk of damaging the brush, always use the original tips. A total of 2 heads included

7.3 Brushing

- Wet the toothbrush head.
- Apply some toothpaste to the wet brush.
- To avoid splashing toothpaste around, insert the brush in your mouth before turning it on.
- Turn the device on by pressing the On/Off button.
- Start by brushing the outside surfaces of the teeth. Press the brush head slightly against the gums and brush the teeth starting from the gums.
- Brush all teeth one after another, starting from the gum.
- Switch to the inner side and repeat these steps.
- Then brush the backs of your molars carefully.
- Once you are finished brushing, turn off the device with the On/Off button. The toothbrush will switch off and the illumination will go off.
- After finishing brushing, wash the brush head under running water and dry.

Caution!

The toothbrush will switch off after 2 minutes.

At first, your gums may bleed slightly as a result of the new type of irritation. This is normal. However, if this condition persists for more than one week, consult your dentist.

8. CLEANING AND CARE

Correct and regular cleaning will improve the safety and extend the life span of the device.

Warning! Turn the device off and leave it to cool down before cleaning to avoid electric shock or burn.

Caution! Do not clean the device with any chemicals, alkali, abrasive or disinfecting agents, as they may damage the surface of the device.

- Clean the toothbrush head thoroughly after brushing by rinsing it under running water.
- For the sake of hygiene, it is recommended that the brush head be changed every 3 months.

9. REPAIR

The device does not contain any user-serviceable parts. Do not attempt to repair the device on your own. Always have a specialist perform repairs.

10. STORAGE AND TRANSPORT

Clean the device as described.

Keep the device in a dry and well-ventilated place inaccessible to children.

Protect the device against vibration and shock during transport.

11. DISPOSAL

All the packaging materials are recyclable and are labelled as such. Dispose of the packaging in accordance with local regulations.

Keep the packaging materials out of children's reach, as these materials can pose a hazard.

Correct disposal of the device:

1.  According to the WEEE Directive (2012/19/EU), the crossed-out wheellie bin symbol (shown on the left) is used to label all electric and electronic devices requiring segregation.
2. Do not dispose of the spent product with domestic waste: turn it over to an electric and electronic devices collection and recycling point. The crossed-out wheellie bin symbol placed on the product, instruction manual or package communicates

this requirement.

3. Plastics contained in the device can be recycled in accordance with the specific marking. Thanks to recycling, re-use of materials or other forms of utilizing used equipment, you make a significant contribution to the protection of our natural environment.
4. Information on electric and electronic devices collection locations is available from local government agencies or from the dealer.
5. Spent or fully discharged single-use and rechargeable batteries must be discarded in dedicated labelled bins, turned over to hazardous waste collectors or returned to electric equipment dealers.

12. CE DECLARATION

This device has been designed, manufactured and marketed in compliance with the requirements of the Electromagnetic Compatibility Directive; therefore, the product has been marked with the CE mark and a declaration of conformity has been issued for it, which is made available to market regulators.

13. WARRANTY

In order to submit a complaint about the product, **submit the device to the Customer Service Point at any Biedronka store.**

Should you have any questions or issues related to product operation, send them to the following e-mail address:

- infolinia@vershold.com
- or contact us by phone: +48 667 090 903

Your opinion is important to us. Evaluate our product at:

www.vershold.com/opinie

1. The Warrantor of this product grants a 24-month warranty valid from the date of purchase. If any defect is found, a complaint about the product should be lodged at the point of purchase.
2. A defective/damaged product means a product failing to provide the features described in the User's Manual due to intrinsic properties of the device.
3. The Beneficiary under the warranty shall deliver the device to the Customer Service Point at a given store in order to submit a complaint about the equipment. When asserting your rights under warranty, you need to present proof of purchase (receipt, invoice) together with a description of the device defect in as much detail as possible, in particular of the external manifestations of the defect.
4. The warranty excludes products with physical damage other than that caused by the manufacturer or distributor, including specifically any damage resulting from incorrect usage or force majeure.
5. The warranty does not cover products from which a tamper-proof seal has been removed and products to which repairs, alterations or structural modifications have been made or attempted.
6. The warranty excludes components subject to natural wear and tear during operation.
7. It is recommended that the complete product be submitted together with the complaint so as to facilitate service operations.
8. This warranty for the sold consumer product does not exclude, restrict or suspend any rights of the Buyer under warranty for

physical and legal defects of goods sold pursuant to the Civil Code of 23 April 1964 (Journal of Laws Dz.U.2014.121).

9. The warranty shall be applicable within the territory of the Republic of Poland.

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